



Engraver



Drybox-2



MD-Concert and MD-Concertino




User-friendly accessories specially designed for materialographic specimens

Marking

Marking of specimens for subsequent identification

If a specimen is not properly marked, the preparation efforts may prove a waste of time if the sample can not be properly identified later.

Struers Engraver has a universal engraving point, suited for marking of practically all metals, plastics, ceramics and glass. A tungsten carbide point is available for extra hard material. The Engraver is very simple to operate: You hold the engraver like you would hold a pencil. The vibration amplitude is regulated by means of a control button.

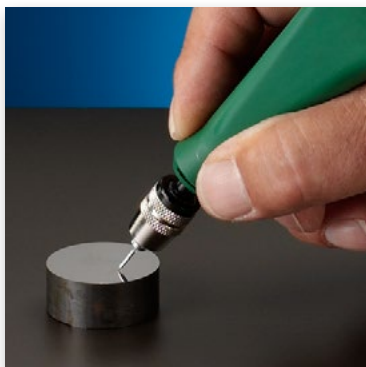

Drying

Drying of specimen holders, single specimens and curing of cold mounting mounts

Struers Drybox-2 is a combined drying and curing apparatus which produces a strong air flow that can be heated between 40 °C and 90°C.

Drybox-2 may considerably reduce the curing time for cold mounted specimens. Epoxy cold mounting systems start to cure due to a chemical reaction when the two components are mixed. Heat advances the start of the reaction, but once initiated, the reaction requires cooling especially for large mounts. Drybox-2 keeps the temperature in the mount at the right level throughout the entire curing process.



Table unit

Storing

Practical, convenient and space saving storage for materialographic laboratories

Storage cabinet for MD-Consumables

The storage cabinets have shelves for storing each individual MD-Consumables, allowing the disc to dry effectively, and protects the preparation surface from contamination. Identification labels to be mounted at the front of each shelf allows individual disc information.

The cabinet version for MD-discs 350 mm dia. has shelves are pulled out from the front. The shelves of the versions for MD-discs 300, 250 and 200 mm dia. swing out either to the left or the right to provide access no matter where the storage cabinet is placed.

Table unit for cut-off machines

Apart from providing support for the machines, this unit comprises a compartment for the recirculation cooling unit so the unit can be rolled out for cleaning. It also has a drawer for tools and shelves for cut-off wheels. An extension for the table unit can be mounted on both right and left hand side. For easy access to the cut-off machine, the table including the machine, can be moved with a lifting truck after mounting a support bar.



Struers A/S
 Pederstrupvej 84
 DK-2750 Ballerup,
 Denmark
 Phone +45 44 600 800
 Fax +45 44 600 801
 struers@struers.dk
 www.struers.com

AUSTRALIA & NEW ZEALAND
Struers Australia
 27 Mayneview Street
 Milton QLD 4064
 Australia
 Phone +61 7 3512 9600
 Fax +61 7 3369 8200
 info.au@struers.dk

JAPAN
Marumoto Struers K.K.
 Takara 3rd Building
 18-6, Higashi Ueno 1-chome
 Taito-ku, Tokyo 110-0015
 Phone +81 3 5688 2914
 Fax +81 3 5688 2927
 struers@struers.co.jp

BELGIQUE (Wallonie)
Struers S.A.S.
 370, rue du Marché Rollay
 F- 94507 Champigny
 sur Marne Cedex
 Téléphone +33 1 5509 1430
 Télécopie +33 1 5509 1449
 struers@struers.fr

NETHERLANDS
Struers GmbH Nederland
 Elektraweg 5
 3144 CB Maassluis
 Telefoon +31 (10) 599 7209
 Fax +31 (10) 5997201
 netherlands@struers.de

BELGIUM (Flanders)
Struers GmbH Nederland
 Elektraweg 5
 3144 CB Maassluis
 Telefoon +31 (10) 599 7209
 Fax +31 (10) 5997201
 netherlands@struers.de

ÖSTERREICH
Struers GmbH
 Zweigniederlassung Österreich
 Betriebsgebiet Puch Nord 8
 5412 Puch
 Telefon +43 6245 70567
 Fax +43 6245 70567-78
 austria@struers.de

CANADA
Struers Ltd.
 7275 West Credit Avenue
 Mississauga, Ontario L5N 5M9
 Phone +1 905-814-8855
 Fax +1 905-814-1440
 info@struers.com

POLAND
Struers Sp. z o.o.
 Oddział w Polsce
 ul. Jasnogórska 44
 31-358 Kraków
 Phone +48 12 661 20 60
 Fax +48 12 626 01 46
 poland@struers.de

CHINA
Struers Ltd.
 No. 1696 Zhang Heng Road
 Zhang Jiang Hi-Tech Park
 Shanghai 201203, P.R. China
 Phone +86 (21) 6035 3900
 Fax +86 (21) 6035 3999
 struers@struers.cn

ROMANIA
Struers GmbH
 Sucursala Sibiu
 Str.Scoala de Inot, nr. 18
 RO-550005 Sibiu
 Phone +40 269 244 558
 Fax +40 269 244 559
 romania@struers.de

CZECH REPUBLIC
Struers GmbH
 Organizační složka
 Havičkova 361
 CZ-252 63 Rostoky u Prahy
 Phone +420 233 312 625
 Fax +420 233 312 640
 czechrepublic@struers.de

SCHWEIZ
Struers GmbH
 Zweigniederlassung Schweiz
 Weissenbrunnstraße 41
 CH-8903 Birmensdorf
 Telefon +41 44 777 63 07
 Fax +41 44 777 63 09
 switzerland@struers.de

DEUTSCHLAND
Struers GmbH
 Carl-Friedrich-Benz-Straße 5
 D-47877 Willich
 Telefon +49 (0) 2154 486-0
 Fax +49 (0) 2154 486-222
 verkauf@struers.de

SINGAPORE
Struers Singapore
 627A Aljunied Road,
 #07-08 BizTech Centre
 Singapore 389842
 Phone +65 6299 2268
 Fax +65 6299 2661
 struers.sg@struers.dk

FRANCE
Struers S.A.S.
 370, rue du Marché Rollay
 F-94507 Champigny
 sur Marne Cedex
 Téléphone +33 1 5509 1430
 Télécopie +33 1 5509 1449
 struers@struers.fr

SUOMI
Struers Suomi
 Kalevankatu 43
 00180 Helsinki
 Puhelin +358 (0)207 919 430
 Faksi +358 (0)207 919 431
 finland@struers.fi

HUNGARY
Struers GmbH
 Magyarországi Fióktelep
 Tatai ut 53
 2821 Gyermely
 Phone +36 (34) 880546
 Fax +36 (34) 880547
 hungary@struers.de

SWEDEN
Struers Sverige
 Ekbacksvägen 22
 168 69 Bromma
 Telefon +46 (0)8 447 53 90
 Telefax +46 (0)8 447 53 99
 info@struers.se

IRELAND
Struers Ltd.
 Unit 11 Evolution@ AMP
 Whittle Way, Catcliffe
 Rotherham S60 5BL
 Tel. +44 0845 604 6664
 Fax +44 0845 604 6651
 info@struers.co.uk

UNITED KINGDOM
Struers Ltd.
 Unit 11 Evolution @ AMP
 Whittle Way, Catcliffe
 Rotherham S60 5BL
 Tel. +44 0845 604 6664
 Fax +44 0845 604 6651
 info@struers.co.uk

ITALY
Struers Italia
 Via Monte Grappa 80/4
 20020 Arese (MI)
 Tel. +39-02/38236281
 Fax +39-02/38236274
 struers.it@struers.it

USA
Struers Inc.
 24766 Detroit Road
 Westlake, OH 44145-1598
 Phone +1 440 871 0071
 Fax +1 440 871 8188
 info@struers.com

Technical Data

Subject	Specification
Struers Engraver	
Strokes per second	100 at 50 Hz power supply 120 at 60 Hz power supply
Drybox-2	
Output	2000 W
Rate of airflow	100 m ³ /hr
Dimensions: (WxDxH)	290 mm x 292 mm x 172 mm
Weight:	4.5 kg
Table Unit	
Dimensions: (WxDxH)	930 mm x 950 mm x 800 mm 36.6" x 37.4" x 31.5"
Extension for table unit	
Dimensions: (WxD)	400 mm x 950 mm 15.7" x 37.4"
Storage Cabinet	
Dimensions: (WxDxH)	379 mm x 397 mm x 240 mm 14.9" x 15.6" x 9.4"

Specifications

	Cat. no.
Struers Engraver For marking of specimens 1 x 110 V / 50-60 Hz 1 x 220 V / 50-60 Hz	00146121 00146133
Engraver Tip Tungsten carbide tipped point for Engraver	40900043
Drybox-2 Combined specimen drier and curing oven with adjustable temperature (approx. 40-90 °C) and low noise fan. With basket for specimens (04236904) and slot for specimen holders 1 x 220 V / 50 Hz 1 x 220 V / 60 Hz	03706133 03706134
Table Unit For Discotom and Labotom cut-off machines. With compartment for recirculation cooling unit, drawer for tools and shelves for cut-off wheels.	06266101
Extension for table unit Can be mounted on both right and left hand side	06266901
Storage Cabinet For 10 MD preparation surfaces in 350 mm (14") dia.	05666001
MD-Concert Storage cabinet for 9 MD preparation surfaces in 300 mm (12"), 250 mm (10") or 200 mm (8") dia.	05306101
MD-Concertino Storage cabinet for 7 MD preparation surfaces in 200 mm (8") dia.	05306102

Struers' products are subject to constant product development. Therefore, we reserve ourselves the right to introduce changes in our products without notice.